

prepFAST CARBON

Status Indicator Light

Alerts user to system problems
(leaks, empty reagents, etc...)

Monitor

Viewing prepFAST CARBON software
with dynamic library

Fluorocapper

Automated bottle cap removal at time of
sampling for 125 mL or 250 mL bottles

Barcode Reading

Identify sample at time of analysis
with barcode reading

Enclosed Ultrapure Fluoronic Autosampler

True fluoropolymer mechanical action
with magnetically coupled linear drive

Autocalibration & Sample Preparation

Autocalibration
Automated MSA
Autodilution

High Purity PFA Reservoirs

Ensures purity of reagents

Integrated Level Sensors

Alerts user when reagents are low

Leak Sensors

Detects leaks and alerts user

Ultra-clean Exhausted Enclosure with Air Shower

Ensures ultra-pure environment
Optional upgrade – ULPA filtered environment

Agilent 6530 QTOF-MS

Infusion mode and speciation mode analyzer
(Instrument included)

HDPE Countertop

Inert and acid resistant

Storage Drawer

Store accessories and consumables

Speciation Mode

Optional add-on for automated
chromatography

Vented Storage Areas


Store bottles, reagents, etc.

Chamber for Vacuum Pump


Pull out drawer for easy maintenance

Demountable Bench

Simplified Installation

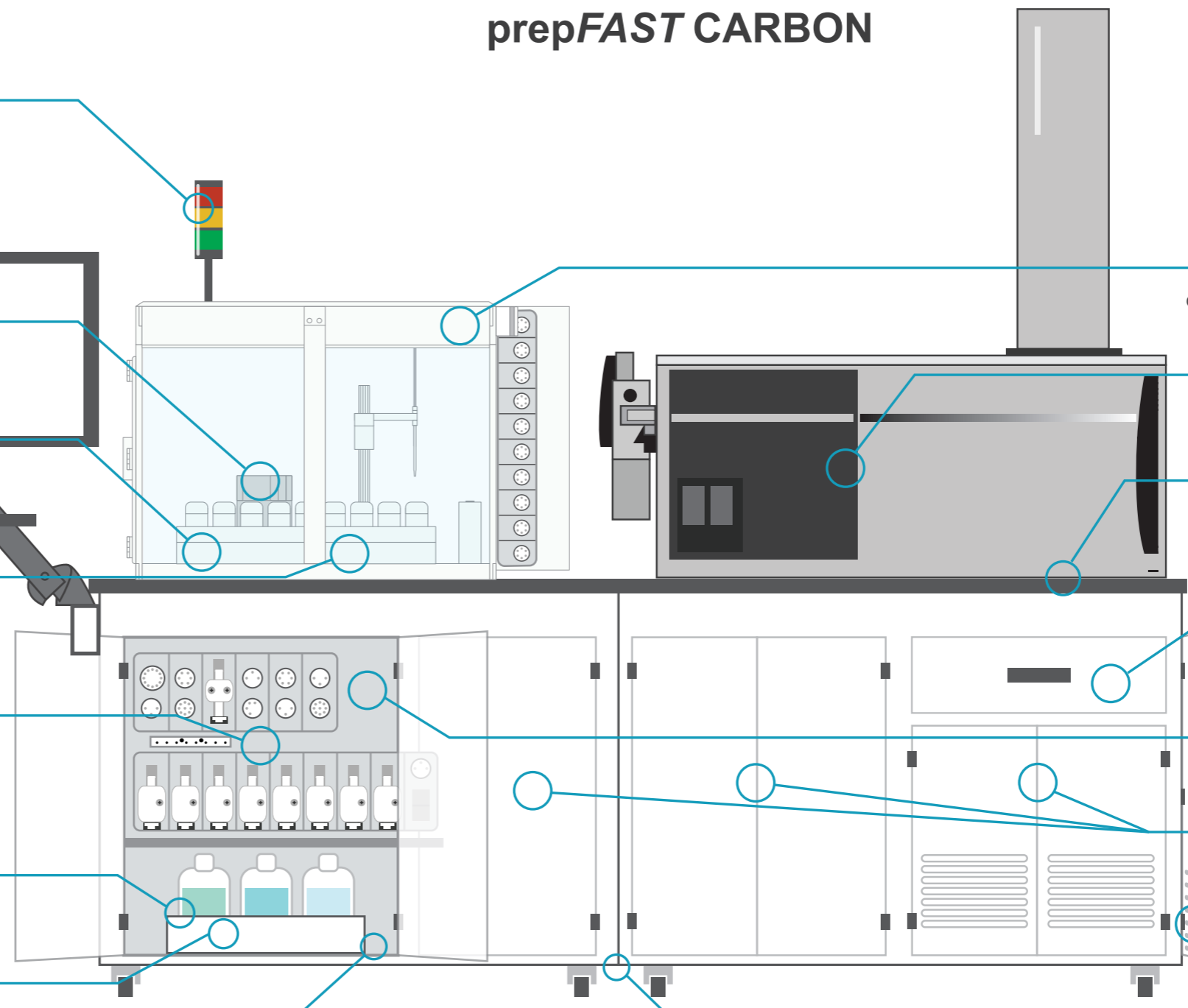

prepFAST CARBON

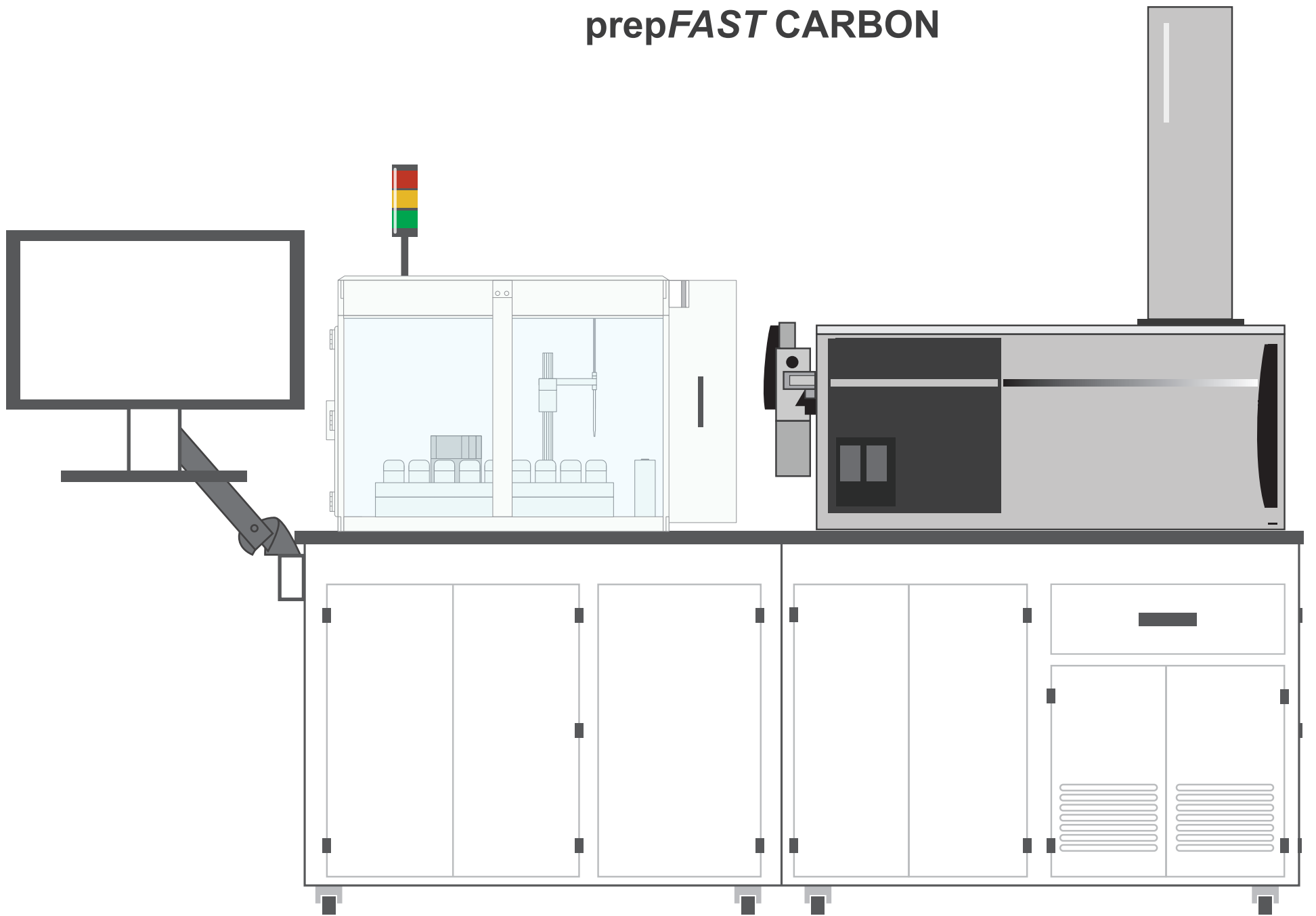

prepFAST CARBON

