

Air Shower

- Customer supplied gas flow provides protective air barrier across openings
- ULPA filtration optional

Ultraclean Enclosure

- Three side easy access doors
- Longitudinal exhaust manifold in sample deck
- Fumes removed at sample level

AutoAlign Arm

- Self-realigns after encountering capped sample tube
- Metal-free fluoropolymer probe aligner

UPW and Calibration Solution

PFA Sample Containers


- <1 mL up to 500 mL

FluroWand Racks

- Custom racks for FluroWand compatible vials, bottles and jars
- Elevated to locate vessel tops at consistent height for FluroWand operation

Syringe Driven

- Precise ultrapure sample uptake and delivery

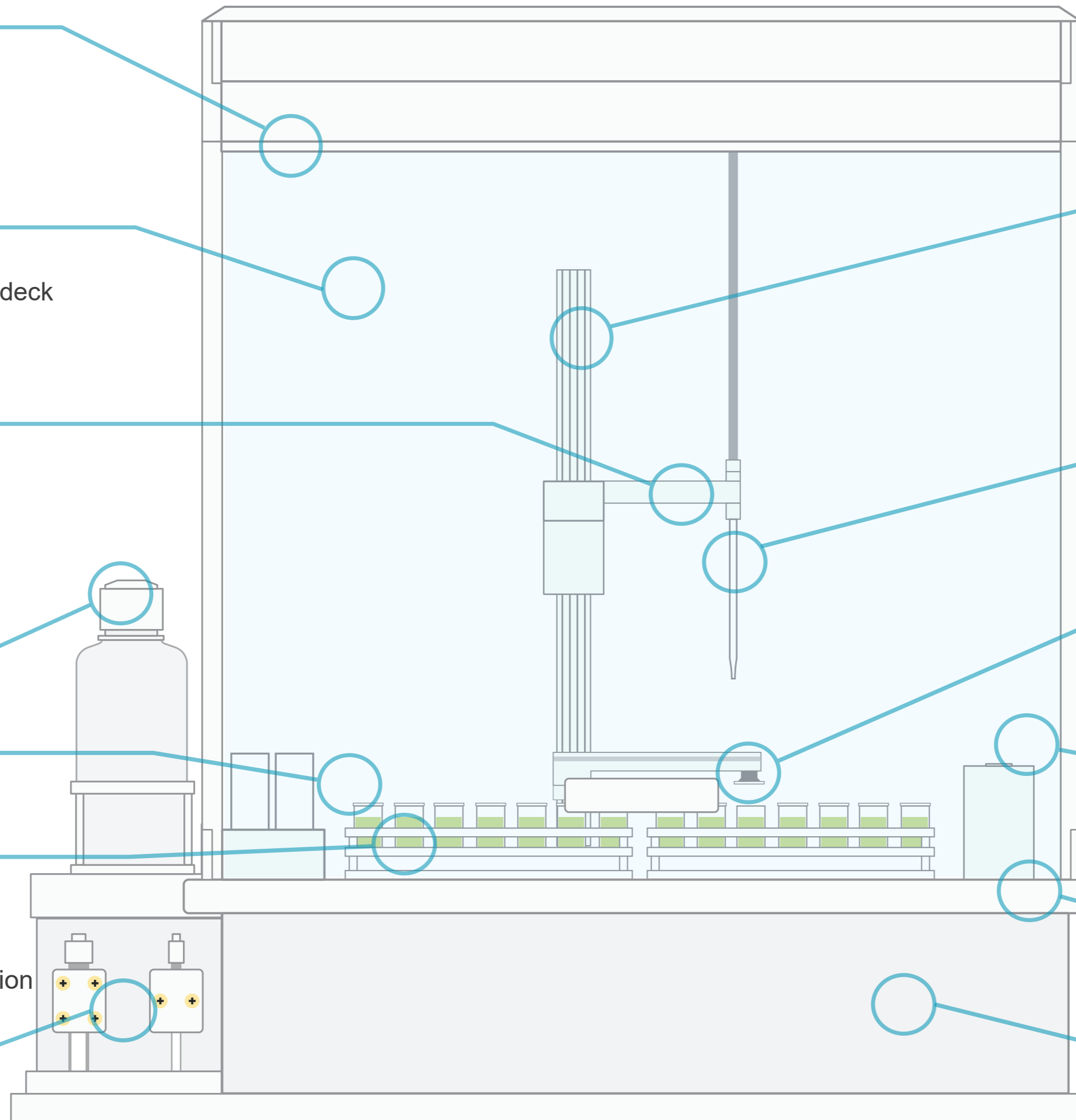

Magnetic Snap Valves

Fluoronic Drive

- True fluoropolymer mechanical action
- Magnetically coupled linear drive
- Chemically resistant & long-lived

FluroProbe

FluroWand Arm

- Removes and replaces sample covers with a vacuum suction cup
- Minimizes risk of contamination
- Prevents sample evaporation and gas phase reactions of samples

Dual Pumped Rinse Station

- Ultra-low cross-contamination
- Independently pumped for low rinse consumption

PTFE Autosampler Deck

Optional Flowing UPW